

IAEA Analytical Quality in Nuclear Applications Series No. 14

Reference Material IAEA 413: Major, Minor and Trace Elements in Algae



IAEA

International Atomic Energy Agency

REFERENCE MATERIAL IAEA 413:
MAJOR, MINOR AND TRACE ELEMENTS IN ALGAE

The following States are Members of the International Atomic Energy Agency:

| | | |
|-------------------------------------|---------------------------|--|
| AFGHANISTAN | GHANA | NORWAY |
| ALBANIA | GREECE | OMAN |
| ALGERIA | GUATEMALA | PAKISTAN |
| ANGOLA | HAITI | PALAU |
| ARGENTINA | HOLY SEE | PANAMA |
| ARMENIA | HONDURAS | PARAGUAY |
| AUSTRALIA | HUNGARY | PERU |
| AUSTRIA | ICELAND | PHILIPPINES |
| AZERBAIJAN | INDIA | POLAND |
| BAHRAIN | INDONESIA | PORTUGAL |
| BANGLADESH | IRAN, ISLAMIC REPUBLIC OF | QATAR |
| BELARUS | IRAQ | REPUBLIC OF MOLDOVA |
| BELGIUM | IRELAND | ROMANIA |
| BELIZE | ISRAEL | RUSSIAN FEDERATION |
| BENIN | ITALY | SAUDI ARABIA |
| BOLIVIA | JAMAICA | SENEGAL |
| BOSNIA AND HERZEGOVINA | JAPAN | SERBIA |
| BOTSWANA | JORDAN | SEYCHELLES |
| BRAZIL | KAZAKHSTAN | SIERRA LEONE |
| BULGARIA | KENYA | SINGAPORE |
| BURKINA FASO | KOREA, REPUBLIC OF | SLOVAKIA |
| BURUNDI | KUWAIT | SLOVENIA |
| CAMBODIA | KYRGYZSTAN | SOUTH AFRICA |
| CAMEROON | LATVIA | SPAIN |
| CANADA | LEBANON | SRI LANKA |
| CENTRAL AFRICAN REPUBLIC | LESOTHO | SUDAN |
| CHAD | LIBERIA | SWEDEN |
| CHILE | LIBYAN ARAB JAMAHIRIYA | SWITZERLAND |
| CHINA | LIECHTENSTEIN | SYRIAN ARAB REPUBLIC |
| COLOMBIA | LITHUANIA | TAJIKISTAN |
| CONGO | LUXEMBOURG | THAILAND |
| COSTA RICA | MADAGASCAR | THE FORMER YUGOSLAV REPUBLIC OF MACEDONIA |
| CÔTE D'IVOIRE | MALAWI | TUNISIA |
| CROATIA | MALAYSIA | TURKEY |
| CUBA | MALI | UGANDA |
| CYPRUS | MALTA | UKRAINE |
| CZECH REPUBLIC | MARSHALL ISLANDS | UNITED ARAB EMIRATES |
| DEMOCRATIC REPUBLIC OF THE CONGO | MAURITANIA | UNITED KINGDOM OF GREAT BRITAIN AND NORTHERN IRELAND |
| DENMARK | MAURITIUS | UNITED REPUBLIC OF TANZANIA |
| DOMINICAN REPUBLIC | MEXICO | UNITED STATES OF AMERICA |
| ECUADOR | MONACO | URUGUAY |
| EGYPT | MONGOLIA | UZBEKISTAN |
| EL SALVADOR | MONTENEGRO | VENEZUELA |
| ERITREA | MOROCCO | VIETNAM |
| ESTONIA | MOZAMBIQUE | YEMEN |
| ETHIOPIA | MYANMAR | ZAMBIA |
| FINLAND | NAMIBIA | ZIMBABWE |
| FRANCE | NEPAL | |
| GABON | NETHERLANDS | |
| GEORGIA | NEW ZEALAND | |
| GERMANY | NICARAGUA | |
| | NIGER | |
| | NIGERIA | |

The Agency's Statute was approved on 23 October 1956 by the Conference on the Statute of the IAEA held at United Nations Headquarters, New York; it entered into force on 29 July 1957. The Headquarters of the Agency are situated in Vienna. Its principal objective is "to accelerate and enlarge the contribution of atomic energy to peace, health and prosperity throughout the world".

IAEA Analytical Quality in Nuclear Applications No. IAEA/AQ/14

***Reference Material IAEA 413:
Major, Minor and Trace Elements
in Algae***

COPYRIGHT NOTICE

All IAEA scientific and technical publications are protected by the terms of the Universal Copyright Convention as adopted in 1952 (Berne) and as revised in 1972 (Paris). The copyright has since been extended by the World Intellectual Property Organization (Geneva) to include electronic and virtual intellectual property. Permission to use whole or parts of texts contained in IAEA publications in printed or electronic form must be obtained and is usually subject to royalty agreements. Proposals for non-commercial reproductions and translations are welcomed and considered on a case-by-case basis. Enquiries should be addressed to the IAEA Publishing Section at:

Sales and Promotion, Publishing Section
International Atomic Energy Agency
Vienna International Centre
PO Box 100
1400 Vienna, Austria
fax: +43 1 2600 29302
tel.: +43 1 2600 22417
email: sales.publications@iaea.org
<http://www.iaea.org/books>

For further information on this publication, please contact:

Terrestrial Environment Laboratory, Seibersdorf
International Atomic Energy Agency
2444 Seibersdorf
Austria

REFERENCE MATERIAL IAEA 413:
MAJOR, MINOR AND TRACE ELEMENTS IN ALGAE
IAEA, VIENNA, 2010
IAEA/AQ/14
ISSN 2074-7659
© IAEA, 2010
Printed by the IAEA in Austria
March 2010

FOREWORD

Reference materials are a basic requirement for any sort of quantitative chemical and radiochemical analysis. Laboratories need them for calibration and quality control throughout their analytical work.

The IAEA started to produce reference materials in the early 1960s to meet the needs of the analytical laboratories in its Member States that required reference materials for quality control of their measurements. The initial efforts were focused on the preparation of environmental reference materials containing anthropogenic radionuclides for use by those laboratories employing nuclear analytical techniques. These reference materials were characterized for their radionuclide content through interlaboratory comparison involving a core group of some 10 to 20 specialist laboratories. The success of these early exercises led the IAEA to extend its activities to encompass both terrestrial and marine reference materials containing primordial radionuclides and trace elements.

Today, the IAEA has more than 90 reference materials and maintains a customer base of about 5000 members from more than 85 Member States.

Within the frame of IAEA activities in production and certification of RM, this report describes the certification of the IAEA 413: Major, minor and trace elements in algae. Details are given on methodologies and data evaluation.

The IAEA officer responsible for this publication was A. Shakhashiro of the IAEA Environment Laboratories.

EDITORIAL NOTE

The use of particular designations of countries or territories does not imply any judgement by the publisher, the IAEA, as to the legal status of such countries or territories, of their authorities and institutions or of the delimitation of their boundaries.

The mention of names of specific companies or products (whether or not indicated as registered) does not imply any intention to infringe proprietary rights, nor should it be construed as an endorsement or recommendation on the part of the IAEA.

CONTENTS

| | | |
|------|--|----|
| 1. | INTRODUCTION | 1 |
| 2. | DESCRIPTION OF ALGAE MATERIAL | 2 |
| 2.1. | Homogeneity study | 2 |
| 2.2. | Long-term stability | 4 |
| 3. | CHARACTERIZATION | 6 |
| 4. | EVALUATION OF RESULTS | 13 |
| 4.1. | General | 13 |
| 4.2. | Statistical screening of the combined data sets | 13 |
| 4.3. | Calculation of property value and associated uncertainty | 15 |
| 4.4. | Traceability of results | 17 |
| 4.5. | Intended use | 17 |
| 4.6. | Instructions for use | 17 |
| 4.7. | Dry mass determination | 18 |
| | APPENDIX: GRAPHICAL PRESENTATION OF THE CHARACTERIZATION RESULTS OF DIFFERENT ANALYTICAL TECHNIQUES | 19 |
| | REFERENCES | 27 |
| | ABBREVIATIONS | 29 |
| | PARTICIPATING LABORATORIES | 31 |
| | CONTRIBUTORS TO DRAFTING AND REVIEW | 33 |

1. INTRODUCTION

A reliable determination of toxic elements in environmental samples plays a vital role in the assessment of environmental pollution. Appropriate assessment depends on the careful selection of representative samples that truly reflect pollution levels. Biomonitors, such as lichen and algae, are examples of environmental samples that have been widely used by the scientific community to monitor environmental pollution. The second pillar supporting reliability of environmental measurements is the availability of appropriate reference materials for analytical quality control and method validation. Besides IAEA-392, IAEA 413 is another algal material prepared for this aim.

Reference values for the mass fractions and associated standard uncertainties have been established for: As (127 ± 6.6 mg/kg), Ca (3143 ± 112 mg/kg), Cd (204 ± 8.5 mg/kg), Co (4.24 ± 0.25 mg/kg), Cr (377 ± 14 mg/kg), Fe (1370 ± 39 mg/kg), K (10740 ± 270 mg/kg), Mg (4058 ± 117 mg/kg), Mn (158 ± 3.4 mg/kg), Na (375 ± 20 mg/kg), Ni (113 ± 4.9 mg/kg), Pb (242 ± 7 mg/kg), Zn (169 ± 3.3 mg/kg).

For Cu (11.1 ± 1.0 mg/kg) and Hg (53.2 ± 4.0 mg/kg) information values are reported.

During sample production and certification, the requirements for reference material production and certification as stated in ISO guides 34 and 35 [1, 2] were taken into account. This report summarizes the preparation and certification process.

The input, support and fruitful discussion with all analysts: E. de Nadai Fernandes, P. Ostapczuk, E. Benetka, P. Bode, M. Campbell, S. Bamfort and M. Makarewicz, as well as with U. Sansone, and P. Martin is highly appreciated.

2. DESCRIPTION OF ALGAE MATERIAL

The IAEA-413 algae material was prepared by the Terrestrial Environment Laboratory of the International Atomic Energy Agency's in cooperation with the Institute of Microbiology, Academy of Sciences of the Czech Republic in Trebon.

The IAEA-413 algae material (type: *Chlorella Boehm*) was produced under standard outdoor culture conditions. The sloped bioreactor bed was constructed from steel. Mineral nutrients were prepared from chemicals and added as required for an optimal growing. The amount of elements added to the nutrient solution was adjusted to allowing bioaccumulation without poisoning the algae [3, 4].

After harvesting, the algae were transferred to a cooled tank and continuously spray dried with moderate venting and at a moderate temperature, avoiding damage to the cell structure.

A two kilograms sample of the dry bulk algae material was irradiated with γ rays to improve long-term stability of the material by reducing microbiological action. Based on the particle size measurements it could be verified that the single cell structure of the algae had been maintained during preparation. Although the natural single cell particle size was around 4.6 μm , the average particle size measured using the Martin's diameter method increased to about 8 or 16 μm due to conglomeration.

Bottling of IAEA-413 was done under normal laboratory conditions, taking all precautions to avoid segregation. After bottling, the samples were packed into plastic bags, sealed, and sterilized by gamma ray irradiation with a total dose of 25 kGy using a Co-60 source.

2.1. Homogeneity study

The homogeneity of the material was tested twice, directly after bottling and during the characterization exercise. Sample bottles covering the whole bottling range were randomly selected for homogeneity tests and the characterization. During the homogeneity tests emphasis was placed on performing measurements under good repeatability conditions.

Table 1 shows the relative repeatability standard deviation (RSD%) of INAA measurements obtained for 15 bottles. The bottle selection was random and covered the whole bottle number range. The sample mass used as test portion was 10 mg.

TABLE 1: RELATIVE STANDARD DEVIATIONS OBTAINED DURING IAEA-413 HOMOGENEITY MEASUREMENTS

| Element | As | Br | Cd | Co | Cr | Fe | Hg | Na | Zn |
|---------|------|------|------|------|------|------|------|------|------|
| RSD %: | 1.17 | 1.18 | 1.24 | 2.49 | 0.94 | 2.25 | 1.55 | 1.13 | 1.74 |

The second homogeneity tests for IAEA 413 covered all elements for which reference and information values have been assigned. The results reported in Table 2 were derived from 2 to 3 sub-samples from 6 bottles analysed under repeatability conditions during the characterization campaign. The sample portion was in all cases about 100 mg or less.

The standard uncertainty associated with the material heterogeneity was calculated using the formulas stated in ISO Guide 35 [2]. Two approaches, ANOVA (formulas 1 and 2), and ‘taking into account insufficiently good measurement repeatability’ (formula 3) were used. The latter calculation method was used in addition since some bottles gave negative values for the between-bottle variance.

$$s_{wb}^2 = MS_{within} \quad (1)$$

$$s_{bb}^2 = \frac{MS_{among} - MS_{within}}{n_0} \quad (2)$$

$$u_{bb}^* = \sqrt{\frac{MS_{within}}{n}} \sqrt[4]{\frac{2}{v_{MS_{within}}}} \quad (3)$$

The abbreviations are explained in the list of abbreviations.

The uncertainties due to sample heterogeneity were generally very small and confirmed the results from the first INAA data (Table1).

TABLE 2: GENERAL STATISTICS AND THE UNCERTAINTY DUE TO HETEROGENEITY TAKING INTO ACCOUNT A SINGLE METHOD AND ABOUT 100 mg SAMPLE MASS

| Element | n_b/n_R | s.e.m. [%] | F_{test} | Method applied | s_{bb} [%] | u_{bb}^* [%] | Reported method u[%] k=1 |
|---------|-----------|------------|------------|----------------|--------------|----------------|--------------------------|
| As | 6/3 | 0.25 | yes (99%) | INAA | 0.55 | 0.39 | 1.64 |
| Ca | 6/2 | 0.21 | yes | INAA | neg | 1.42 | 1.88 |
| Cd | 6/3 | 0.10 | yes | INAA | 0.80 | 2.86 | 2.40 |
| Co | 6/3 | 0.07 | yes | INAA | 0.30 | 0.46 | 1.20 |
| Cr | 6/3 | 0.05 | yes | INAA | 0.55 | 0.30 | 1.45 |
| Cu | 6/2 | 0.28 | yes | ICP-MS | 2.18 | 1.42 | no info |
| Fe | 6/3 | 0.09 | yes | INAA | 0.83 | 0.53 | 1.34 |
| Hg | 6/3 | 0.15 | yes | INAA | 1.62 | 0.79 | 2.02 |
| K | 6/2 | 0.14 | yes | INAA | 0.51 | 0.87 | 0.94 |
| Mg | 6/2 | 0.06 | yes | ICP-MS | 0.28 | 0.37 | no info |

| Element | n_b/n_R | s.e.m. [%] | F_{test} | Method applied | s_{bb} [%] | u^*_{bb} [%] | Reported method u [%] $k=1$ |
|---------|-----------|------------|-------------------|----------------|--------------|----------------|-------------------------------|
| Mn | 6/2 | 0.11 | yes | INAA | 0.55 | 0.62 | 0.65 |
| Na | 6/3 | 0.05 | yes | INAA | neg | 0.34 | 2.37 |
| Ni | 6/3 | 0.14 | yes | INAA | 0.80 | 0.89 | 2.06 |
| Pb | 6/2 | 0.09 | yes | ICP-MS | 0.64 | 0.48 | 1.31 |
| Zn | 6/3 | 0.09 | yes | INAA | 0.80 | 0.49 | 1.32 |

The facts that the standard errors of the mean (s.e.m.) were below 0.28 % for all elements, and all data sets of the homogeneity study showed normal distribution and had no outlying results, gave additional evidence on the good homogeneity of the material. These data also confirmed that the analytical methods were under control and suitable for the homogeneity test.

The uncertainties associated with the sample heterogeneity (s_{bb} and u^*_{bb}) reported in Table 2 can be used as additional information. They are referring to a defined test portion measured with a stated analysis technique while the certification is taking into account the in-bottle and between-bottle variances of all participating laboratories (see section 5: Evaluation of results). The component of the combined uncertainty of the certification associated with the material heterogeneity cannot be linked to a specific mass of the test portion but is covering a mass range.

During a co-ordinated research project (CRP) on “Reference Materials for Microanalytical Nuclear Techniques” it could be demonstrated that IAEA 413 Algae showed good homogeneity for K, Cr, Mn, Fe, Zn, As, Br and Cd even for a 1 mg test portion [5,6].

Although there are indications that due to the natural small particle size of this material the homogeneity of the material is generally suitable even for smaller test portions, it is recommended that the minimum test portion used for analysis should be 100 mg. If smaller test portion is taken the uncertainty of the property value should be expanded taking into account the relationship of mass and heterogeneity explained in the concept of Ingamels sampling constant or other related concepts [7,8,9].

2.2. Long-term stability

The long-term stability was measured for 10 elements. Storage temperatures of 40, 20, 4 and -20 °C were tested. Test time was about 3 years and ten ($n = 10$) independent determinations were performed for the four temperature conditions.

Applying the F-test and t-tests for the 10 results, the long term reproducibility showed equal variances for all four data sets. In addition the t-tests assured that there was no significant difference on the analysis results of all tested elements for the four storage conditions.

For the trend analysis, the slope (b_1), intercept (b_0) and the uncertainties of both parameters were calculated (s_{b_1} , s_{b_0}). Table 3 shows the results of the trend analysis of samples stored at room

temperature (about 20 °C ± 5 °C) during 1027 days. Five measurements were performed during the first 100 days, and then the measurement intervals were extended to approximately 6 months.

For the calculation of all parameters the equations as described in ISO Guide 35 [2] Chapter 8.3.1 were used. The critical value was calculated:

$$\text{Critical value} = t_{0.95, n-2} * sb_1 \quad (4)$$

The trend would be significant if the slope (b_1) was bigger than the critical value. Since the critical value was always higher, the trend is insignificant. As a consequence it may be stated that no degradation was observed.

TABLE 3: RESULTS OF THE TREND CALCULATION FOR IAEA 413

| | n | b_1 | sb_1 | b_0 | sb_0 | Critical value (95% conf.) | Significant trend |
|----|----|---------|--------|--------|--------|-------------------------------|-------------------|
| Fe | 10 | 0.0293 | 0.1159 | 1453.2 | 48.17 | 0.2674 | no |
| Cu | 10 | 0.0006 | 0.0009 | 10.70 | 0.39 | 0.0020 | no |
| Co | 10 | 0.0002 | 0.0004 | 4.16 | 0.15 | 0.0008 | no |
| Cr | 10 | 0.0191 | 0.0184 | 372.1 | 7.6 | 0.0423 | no |
| Mn | 10 | -0.0051 | 0.0125 | 167.7 | 5.2 | 0.0289 | no |
| Ni | 10 | 0.0003 | 0.0050 | 114.7 | 2.1 | 0.0115 | no |
| Zn | 10 | 0.0086 | 0.0139 | 163.6 | 5.8 | 0.0320 | no |
| As | 10 | 0.0047 | 0.0060 | 126.3 | 2.5 | 0.0139 | no |
| Cd | 10 | -0.0083 | 0.0111 | 203.8 | 4.6 | 0.0257 | no |
| Pb | 10 | 0.0163 | 0.0094 | 245.1 | 3.1 | 0.0262 | no |

It was decided not to include the uncertainty related to the stability study, since the measurement reproducibility standard deviation within the three years of the study was large and would lead to an overestimation of the combined uncertainty. The long-term stability of the material is further monitored using an isochronic approach.

3. CHARACTERIZATION

For the characterization exercise, laboratories were selected based on their successful participation in previous IAEA intercomparison runs. All of them had an implemented quality assurance system conforming to ISO 17025, but not all were formally accredited at the time of the exercise.

During the planning meeting, a short guideline was prepared where quality aspects (e.g. traceability of standards and calibration, instrument and balance performance, quality of pipettes and volumetric flasks) and analytical requirements (e.g. number of sub-samples, test portion mass, dissolution techniques) were defined [10].

Each analyst had received 6 bottles of the algae material (randomly selected covering the whole bottling sequence) and was requested to analyze from each bottle at least 2 sub-samples following the established guidelines.

Test portion mass for the analysis was proposed to be 100 mg (if feasible, to achieve a good limit of detection). To assess whether some elements are silicate bound, the IAEA laboratories performed two dissolution techniques one of them applying HF. Both sets of dissolved samples were measured under repeatability conditions. Results obtained using ICP-OES are shown in Table 4. No significant difference (95% conf.) of the elemental mass fractions measured could be detected.

TABLE 4: ICP-OES RESULTS USING TWO DIFFERENT DIGESTION PROCEDURES (*FOR POTASSIUM INSTEAD OF HClO₄, H₂O₂ AND AAS WAS USED)

| Element | PMW HNO ₃ , HClO ₄ , HF | | | PMW HNO ₃ , HClO ₄ | | | Δ | $\pm\Delta\%$ |
|---------|---|-------|-----|--|-------|-----|----------|---------------|
| | Mass fraction [mg/kg] | \pm | u | Mass fraction [mg/kg] | \pm | u | | |
| As | 126 | \pm | 6 | 127 | \pm | 5 | 0.9 | 0.70 |
| Ca | 3121 | \pm | 87 | 3133 | \pm | 111 | 12 | 0.38 |
| Cd | 200 | \pm | 3 | 201 | \pm | 3 | 1 | 0.52 |
| Co | 4.1 | \pm | 0.3 | 4.0 | \pm | 0.3 | 0.1 | 2.6 |
| Cr | 373 | \pm | 6 | 377 | \pm | 7 | 4 | 1.0 |
| Cu | 10.8 | \pm | 0.4 | 10.7 | \pm | 0.6 | 0.1 | 1.2 |
| Fe | 1359 | \pm | 28 | 1372 | \pm | 22 | 13 | 0.9 |
| Hg | No data available | | | | | | | |
| K* | 10686 | \pm | 171 | 10697 | \pm | 174 | 11 | 0.11 |

| Element | PMW HNO ₃ , HClO ₄ , HF | | | PMW HNO ₃ , HClO ₄ | | | Δ | ±Δ% |
|---------|---|---|----|--|---|-----|---|------|
| | Mass fraction ± u [mg/kg] | | | Mass fraction ± u [mg/kg] | | | | |
| Mg | 4122 | ± | 59 | 4121 | ± | 103 | 1 | 0.03 |
| Mn | 155 | ± | 2 | 157 | ± | 2 | 2 | 0.89 |
| Ni | 110 | ± | 3 | 111 | ± | 2 | 1 | 1.1 |
| Pb | 231 | ± | 5 | 233 | ± | 4 | 2 | 0.56 |
| Zn | 165 | ± | 5 | 165 | ± | 4 | 0 | 0 |

Tables 5.1 to 5.15 give a short summary for all elements on the dissolution procedures used before analysis, and the analytical method applied. The abbreviations are explained on page 25.

It is worth mentioning that all digestion methods were using strong oxidizing acids, acid mixtures or pressure, or a combination of all.

TABLE 5.1 ARSENIC

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|--|-----------------------------------|
| P6/5 | 200 | HPD at 300 °C, HNO ₃ | HG- AAS, 193.7 nm |
| P2/1a | 100 | PMW HNO ₃ , H ₂ O ₂ | ICP-MS, Mass 75 |
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 193.695 nm |
| ND/1 | NA | Non destructive | ED-XRF, As-K _α |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES axial, 188.979 nm |
| ND/2 | 100 | Non destructive | INAA, ⁷⁶ As at 559 keV |
| P2/1 | 100 | PMW HNO ₃ , H ₂ O ₂ | ET AAS , 193.7 nm |
| ND/3 | 150 | Non destructive | INAA ¹ [11] |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, Mass 75 |

¹ Laboratory ND/3 used for the identification and quantification of each element all gamma ray lines, including sum-peaks of all radionuclides produced upon reactor neutron irradiation from the respective stable isotopes of that element. The radionuclides, gamma-ray lines and intensity ratios are traceable to M.Blaauw [11].

TABLE 5.2 CALCIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------------|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 393.367 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 422.673 nm |
| ND/1 | NA | Non destructive | ED-XRF, Ca-K _α |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 422.7 nm |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES axial 422.673 nm |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 44 |
| ND/3 | 150 | Non destructive | INAA [11] |

TABLE 5.3 CADMIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|--|
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 226.502 nm |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 114 |
| ND/3 | 150 | Non destructive | INAA [11] |
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 226.502 nm |
| ND/1 | 100 | Non destructive | INAA, ¹¹⁵ Cd, 336.2 and 527.9 keV |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS |
| ND/2 | 100 | Non destructive | INAA, ¹¹⁵ Cd, 336.2 and 527.9 keV |
| ND/1 | NA | Non destructive | ED-XRF, Cd-K _α |
| P1/1 | 100 | PMW, HNO ₃ | ICP-MS, mass 111 |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 228.802nm |

TABLE 5.4 COBALT

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|--|---|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 228.616 nm |
| P6/5 | 200 | HPD at 300 °C, HNO ₃ | ASV (Square Wave Detection) |
| ND/2 | 100 | Non destructive | INAA, ⁶⁰ Co, 1173 and 1333 keV |
| ND/1 | 100 | Non destructive | INAA, ⁶⁰ Co, 1173 and 1333 keV |
| ND/3 | 150 | Non destructive | INAA* [11] |
| P2/1a | 100 | PMW HNO ₃ , H ₂ O ₂ | ICP-MS, mass 59 |
| P2/1 | 100 | PMW HNO ₃ , H ₂ O ₂ | ET AAS, 242.5nm |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 59 |

TABLE 5.5 CHROMIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|--|-----------------------------------|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 267.716 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 267.716 nm |
| ND/2 | 100 | Non destructive | INAA, ⁵¹ Cr, 320.1 keV |
| ND/1 | 100 | Non destructive | INAA, ⁵¹ Cr, 320.1 keV |
| ND/1 | NA | Non destructive | ED-XRF, Cr-K _α |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 205.552 nm |
| P2/1a | 100 | PMW HNO ₃ , H ₂ O ₂ | ICP-MS, mass 53 |
| P2/1 | 100 | PMW HNO ₃ , H ₂ O ₂ | Flame AAS, 242.5 nm |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 53 |

TABLE 5.6 COPPER

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------------|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 324.754 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 327.396 |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | ET AAS, 324.8 nm |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES axial, 327.396 nm |
| P1HF/1 | 100 | PMW, HNO ₃ , HF | ICP-MS, mass 63 |

TABLE 5.7 IRON

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 259.940 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 259.940 nm |
| ND/2 | 100 | Non destructive | INAA, ⁵⁹ Fe, 1099 and 1292 keV |
| P2/1 | 200 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 248.3 nm |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 57 |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 238.204 nm |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 57 |
| ND/1 | 100 | Non destructive | INAA, ⁵⁹ Fe, 1099 and 1292 keV |

TABLE 5.8 POTASSIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------------------|
| ND/1 | NA | Non destructive | ED-XRF, K-K _α |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 769.898 nm |
| ND/2 | 100 | Non destructive | INAA, ⁴² K, 1525 keV |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 769.9 nm |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 769.896 nm |

TABLE 5.9 MAGNESIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 279.553 nm |
| P1HF/1 | 100 | PMW, HNO ₃ , HF | ICP-MS, mass 24 |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 279.553 nm |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 202.6 nm |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 285.213 nm |

TABLE 5.10 MANGANESE

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------|
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 257.610 nm |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 55 |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 279.8 nm |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 257.610 nm |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 55 |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 257.610 nm |

TABLE 5.11 MERCURY

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|-----------------------------------|------------------------------------|
| P7/4 | 1000 | reflux, HNO ₃ | CV AAS |
| ND/2 | 100 | Non destructive | INAA, ²⁰³ Hg, 279.2 keV |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P5/5 | 300 | LPD, HNO ₃ | CV AAS, with amalgam accumulation |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS mass 202 |

TABLE 5.12 SODIUM

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|----------------------------------|
| P2/1 | 200 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 589.0 nm |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 589.592 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 589.592 nm |
| ND/2 | 100 | Non destructive | INAA, ²⁴ Na, 1369 keV |

TABLE 5.13 NICKEL

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------------------|
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | ET AAS, 232.0 nm |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 60 |
| ND/1 | 100 | Non destructive | INAA, ⁵⁸ Co, 811 keV |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 60 |
| ND/2 | 100 | Non destructive | INAA, ⁵⁸ Co, 811 keV |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 231.604 nm |
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 221.647 nm |
| P6/5 | 200 | HPD at 300 °C, HNO ₃ | ASV |

TABLE 5.14 LEAD

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|---------------------------|
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | ET AAS, 283.3 nm |
| P1/1 | 100 | PMW, HNO ₃ , | ICP-MS, mass 208 |
| ND/1 | NA | Non destructive | ED-XRF, Pb-L _α |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 208 |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES axial, 231.604 nm |
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 220.353 nm |
| P8/5 | 200 | PTFE at 180 °C, HNO ₃ | ICP-OES, 220.353 nm |

TABLE 5.15 ZINC

| Prep/Lab Code | Approx. test portion mass [mg] | Dissolution technique/ Acid(s) | Analysis technique |
|---------------|--------------------------------|---|----------------------------------|
| P2/1 | 100 | PMW, HNO ₃ , H ₂ O ₂ | Flame AAS, 213.9 nm |
| P3/1 | 1000 | PMW, HNO ₃ , HClO ₄ | ICP-OES, 206.200 nm |
| P2/1a | 100 | PMW, HNO ₃ , H ₂ O ₂ | ICP-MS, mass 64 |
| ND/1 | NA | Non destructive | ED-XRF, Zn-K _α |
| ND/2 | 100 | Non destructive | INAA, ⁶⁵ Zn, 1116 keV |
| P4/4 | 1000 | open, HNO ₃ , HClO ₄ | ICP-OES, 213.856 nm |
| ND/3 | 150 | Non destructive | INAA * [11] |
| P1/1 | 100 | PMW, HNO ₃ | ICP-MS, mass 64 |

4. EVALUATION OF RESULTS

4.1. General

All results were collected and used as a basis for discussion. During the characterization meeting, analytical results were rejected based on the following conditions:

- if there were technically justified doubts on the analytical quality, or
- where analysts had used two different digestion methods to evaluate the complete dissolution of the sample and the measurement technique was the same, the analysts were asked to withdraw one set of data. This was to avoid the double counting of methods using dissolution techniques in the evaluation.

For some elements, non-destructive analytical techniques were not available. But since the dissolution methods were independent and very different, and elements were not silica bound (see Table 4), it is assumed that no bias was introduced which could be based on incomplete dissolution.

4.2. Statistical screening of the combined data sets

Since the number of data sets was small for many elements (e.g. $n=5$), their statistical evaluation was expanded to cover all individual results of all methods and analysts and not only the laboratory mean values. The HISTO program was utilized for this purpose. HISTO is a software package developed by the Chemistry Unit of the Terrestrial Environment Laboratory for statistical evaluation of data. Beside general descriptive statistics, the following tests are included in the program:

- Outlier tests (Dixon, Grubbs, Skewness, Kurtosis)
- Directional tests (Skewness, Kurtosis)
- Normality tests (Kolmogorov-Smirnov, Kolmogorov-Smirnov-Lilliefors)

TABLE 6: INFORMATION ON THE STATISTICAL DATA EVALUATION

| Element | n_L | Method used | n_i | $n_{outlier}$ | Normal distribution | s_i [%] | s.e.m. [%] |
|---------|-------|---|-------|---------------|---------------------|-----------|------------|
| As | 10 | INAA (3), HG-AAS, ICP-OES (2), AAS, XRF, ICP-MS (2) | 137 | 0 | yes (99%) | 5.40 | 0.46 |
| Ca | 7 | INAA, AAS, ICP-OES(3), XRF, ICP-MS | 84 | 0 | yes | 3.26 | 0.36 |
| Cd | 10 | INAA(3), ICP-OES(3), AAS, XRF, ICP-MS(2) | 120 | 0 | yes (99%) | 4.01 | 0.37 |
| Co | 8 | INAA(3), ICP-OES, ASV, AAS, ICP-MS(2) | 114 | 0 | yes (99%) | 4.52 | 0.42 |
| Cr | 10 | INAA(3), ICP-OES, AAS, XRF, ICP-MS(2) | 132 | 0 | yes | 3.28 | 0.29 |
| Cu | 5 | ICP-OES(3), AAS, ICP-MS | 60 | 0 | no | 3.87 | 0.50 |
| Fe | 9 | INAA(3), ICP-OES(3), AAS, ICP-MS(2) | 125 | 2 | yes | 3.47 | 0.31 |
| Hg | 5 | INAA(2), AAS CV(2), ICP-MS | 65 | 1 | yes (99%) | 5.73 | 0.71 |
| K | 6 | INAA(2), ICP-OES(2), AAS, XRF | 78 | 0 | yes | 2.35 | 0.27 |
| Mg | 5 | ICP-OES(3), AAS, ICP-MS | 60 | 0 | yes | 2.50 | 0.32 |
| Mn | 7 | INAA(1), ICP-OES(3), ICP-MS (2), AAS | 84 | 1 | yes | 2.79 | 0.30 |
| Na | 5 | INAA(2), ICP-OES(2), AAS | 66 | 0 | yes (99 %) | 4.16 | 0.51 |
| Ni | 9 | ICP-OES(2), AAS, ICP-MS(2), ASV, INAA(3) | 126 | 0 | yes | 4.94 | 0.44 |
| Pb | 7 | ICP-OES(3), AAS, ICP-MS(2), XRF | 84 | 0 | yes | 2.67 | 0.29 |
| Zn | 8 | INAA(2), ICP-OES(2), AAS, XRF, ICP-MS (2) | 102 | 2 | yes | 2.27 | 0.22 |

The standard deviation (s_i) of all individual measurements (n_i) was below 6 % (average 3.6 %) and standard error of the mean (s.e.m) was below 0.71 %. Nearly no outlying data were observed. The directional tests did not always pass the acceptance criteria but still the Kolmogorov-Smirnov and Kolmogorov-Smirnov-Lilliefors normality tests showed normal distributions of the data sets except for Cu.

4.3. Calculation of property value and associated uncertainty

The property values of all elements were established on the basis of a robust approach proposed by David L. Duewer [14] and the Mixture Model Median (MM-median) of the analytical results reported by the expert laboratories was calculated. The MM-median is a direct analogue of the median. It is the location which divides the Mixture Model Probability Density Function (MM-PDF) into two sections of equal area. The MM-median is closely related to the median. It is robust to outliers and also accounts for the reported uncertainty of each measurement result.

To estimate the standard uncertainty associated with the property value the MM-median based Standard Deviation $S(\text{MM-median})$ was calculated from the span of the central 50% of the MM-PDF density function.

Table 7 shows the assigned property values and associated standard uncertainties.

A ‘reference’ property value was assigned if:

- The analyte results passed the normality tests (at least at 99 % confidence level);
- The number of statistical outliers in the data sets was below 5%;
- The maximum standard uncertainty was below the in-house established limits:

| Mass fraction range [mg/kg] | Maximum u (%) |
|-----------------------------|---------------|
| 1-25 | 10 |
| 25-100 | 8 |
| >100 | 6 |

The variances in Table 7 (within-bottle (S_{wb}^2), between-bottle (S_{bb}^2), and between laboratory (S_{bL}^2)) obtained from the fully nested ANOVA. This is one approach listed in the ISO guide 35. All related formulas can be found in Annex A of the guide [2].

TABLE 7: REFERENCE MASS FRACTIONS AND STANDARD UNCERTAINTIES [mg/kg DRY MASS] ESTABLISHED DURING THE CERTIFICATION OF IAEA 413

| Element | Reference value [mg/kg] | | | | Combined standard uncertainty |
|---------|----------------------------|----------|----------|----------|----------------------------------|
| | | u_{wb} | u_{bb} | u_{bL} | u_x |
| As | 127 | 0.12 | 0.17 | 2.06 | 6.6 |
| Ca | 3143 | 5.3 | 5.8 | 33.4 | 112 |
| Cd | 204 | 0.26 | 0.12 | 2.57 | 8.5 |
| Co | 4.24 | 0.05 | 0.00 | 0.06 | 0.25 |
| Cr | 377 | 0.38 | 0.17 | 3.83 | 14 |
| Fe | 1370 | 1.84 | 1.77 | 14.73 | 39 |
| K | 10740 | 17.3 | 7.8 | 87.2 | 270 |
| Mg | 4058 | 4.8 | 8.3 | 41.1 | 117 |
| Mn | 158 | 0.29 | 0.13 | 1.40 | 3.4 |
| Na | 375 | 0.89 | 0.94 | 6.36 | 20 |
| Ni | 113 | 0.42 | 0.00 | 1.10 | 4.9 |
| Pb | 242 | 0.30 | 0.22 | 2.30 | 7 |
| Zn | 169 | 0.24 | 0.35 | 0.62 | 3.3 |

| Element | Information value[mg/kg] | | | | Combined standard uncertainty |
|---------|-----------------------------|----------|----------|----------|----------------------------------|
| | | u_{wb} | u_{bb} | u_{bL} | u_x |
| Cu | 11.1 | 0.03 | 0.01 | 0.17 | 0.5 |
| Hg | 53.2 | 0.14 | 0.13 | 1.45 | 4.0 |

Where:

$$u_{wb} = \frac{S_{wb}^2}{n_i}, \quad u_{bb} = \frac{S_{bb}^2}{n_b \cdot n_L}, \quad u_{bL} = \frac{S_{bL}^2}{n_L}$$

For Hg no reference but only an information value was assigned since the sample dissolution sometimes showed problems which could not be clearly explained. Cu failed to be a reference value, because the data did not show a normal distribution.

4.4. Traceability of results

The property values assigned to the algae reference material are element mass fractions, expressed in the derived SI unit mg/kg. The utmost care was taken regarding the metrological traceability of the property values assigned to this reference material already at the planning phase and during the entire characterization process. Laboratories participating in the characterization campaign have been requested to carefully choose the calibrants and to provide the IAEA with all related information, including certificates. However, the selection of measurement methods and measurement procedures, as well as respective calibrants, was based on the decision of the participating laboratory. A consequence of the use of different calibrants is the fact that the metrological chain(s) for each of the assigned quantity values (combined from a number of results), cannot easily be described. Therefore, the assigned property values — the element mass fractions — although expressed in the derived SI unit, are not intended for calibration purposes, and the reference material as such is not to be used as a calibrant.

4.5. Intended use

This reference material is intended to be used for quality assurance purposes, basically as a quality control material for the measurement of the elemental composition of biological materials, especially biomonitors, for the assessment of a laboratory's analytical work and for the validation of analytical methods.

The estimated standard uncertainty is relatively large due to the consideration of the between-laboratories dispersion in the calculation of the S(MM-median). It is expected that individual laboratories applying a single analytical procedure will produce results with a smaller dispersion. Therefore it is recommended that the users establish their own reproducibility standard deviations to be used as control limits for precision.

4.6. Instructions for use

Before each sub-sampling, the bottle should be shaken thoroughly for 5 min. to re-homogenize the sample. The recommended minimum test portion size is 100 mg.

Analysts are reminded to take appropriate precautions to avoid contamination of the sample and the remaining material in the bottle.

It is recommended to store the material after the opening of the bottle in a refrigerator (at 4 to 8 °C) or dessicator. Exposure to sunlight should be avoided.

Storage at room temperature and even at temperatures up to 40 °C did not show degradation of the originally sealed algae material.

4.7. Dry mass determination

All reference and information values are expressed on a dry mass basis. Therefore the analytical results need to be corrected for the moisture content of the sample at the time of analysis.

It is recommended to dry a separate sample portion of at least 500 mg for 4 hours at 80 °C. If smaller sample test portions are taken, the uncertainty on the dry mass correction factor is increased and should be taken into account for the total uncertainty calculation. Indication of the magnitude of the dry mass uncertainties and their influence on the final results can be found in a related publication [13]. The use of a different drying procedure may lead to biases in the results and should be avoided [13].

APPENDIX
GRAPHICAL PRESENTATION OF THE CHARACTERIZATION
RESULTS OF DIFFERENT ANALYTICAL TECHNIQUES


Laboratory means are reported for all techniques and laboratories used for the certification, they are the mean of all individual analysis performed in six bottles and reported with 95% confidence on their standard deviation.

Reference and information values are reported with expanded uncertainty U ($k = 2$) estimated as a standard deviation of the Mixture Model Median.


The red line represents the derived property value, the blue lines represent the property value $\pm U$ ($k = 2$).

The information on the x-axis gives the laboratory code and the method used for the characterization. Details on the analysis method can be found in section 5 of the report.


IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value $\pm U$ for As

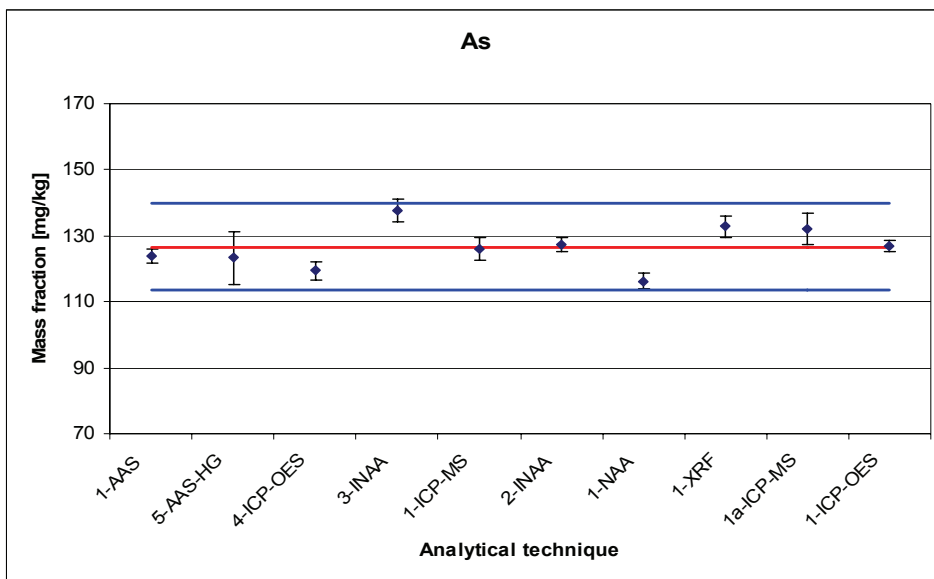


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value $\pm 1\sigma$ for Ca**

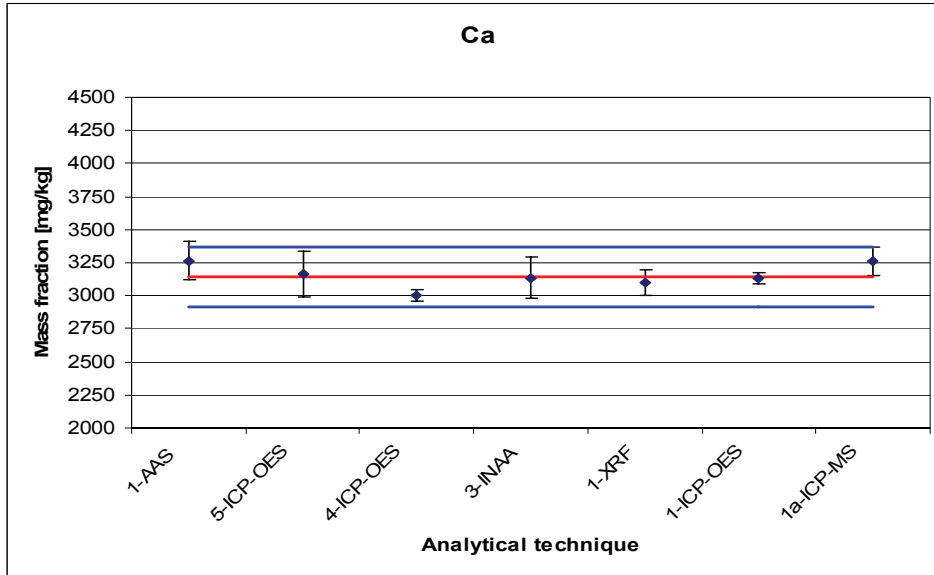


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value $\pm U$ for Cd**

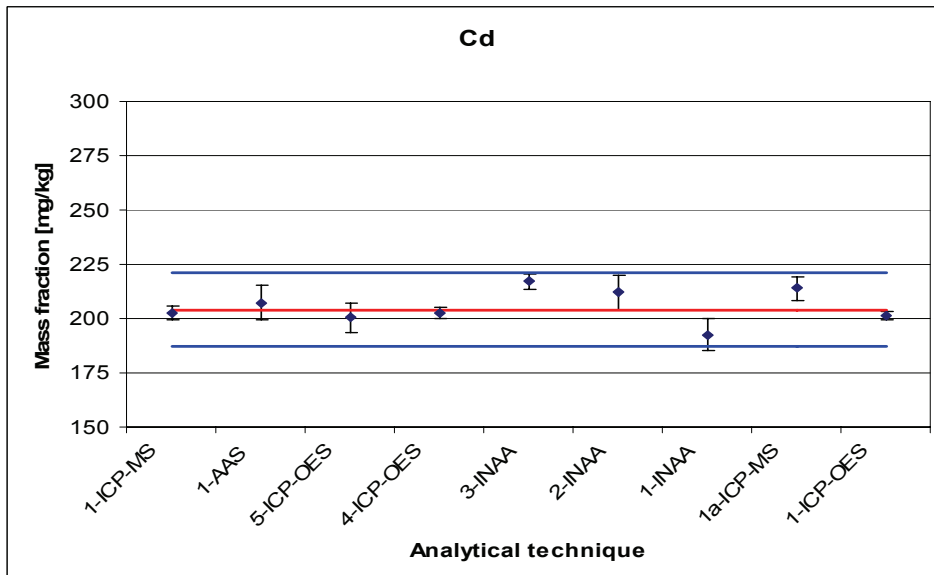


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Co**

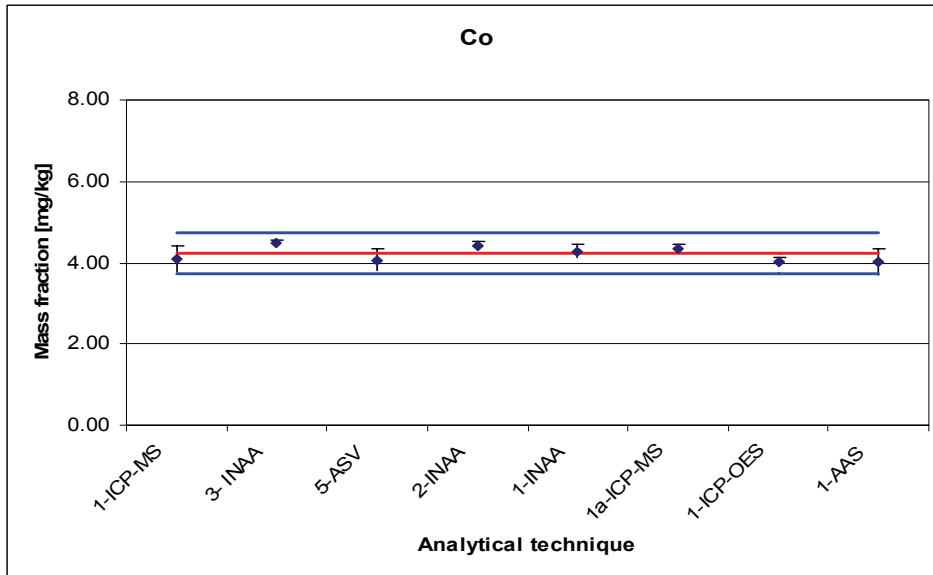


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Cr**

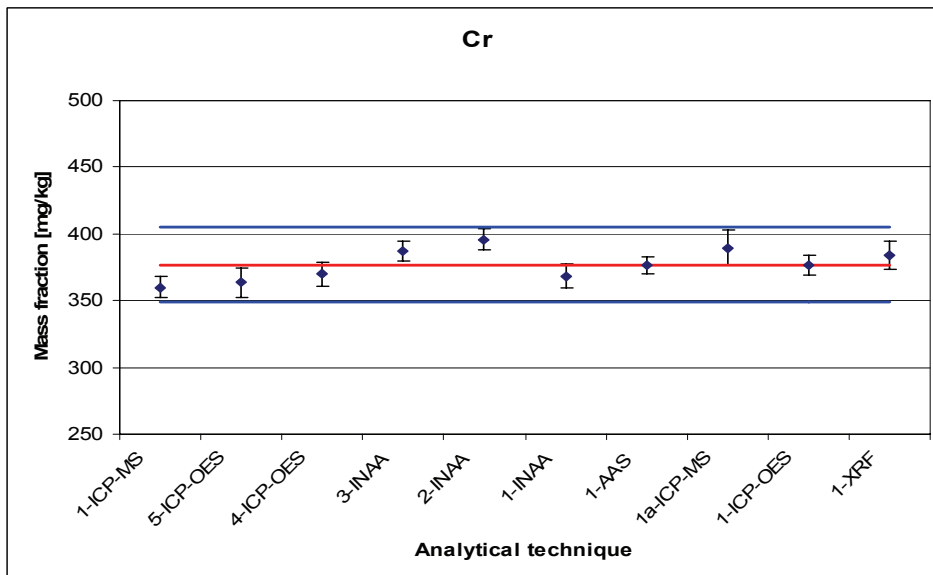


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and information value $\pm U$ for Cu**

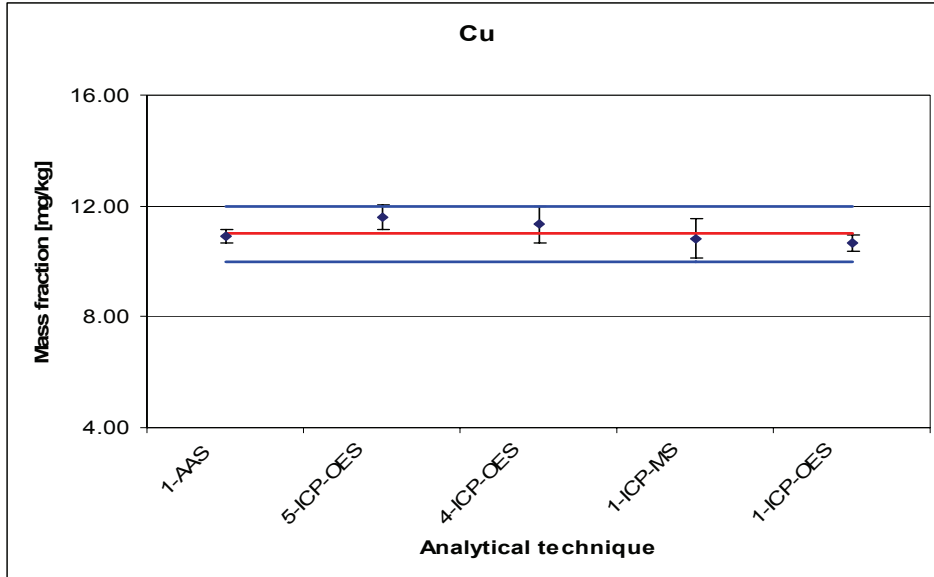


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value $\pm U$ for Fe**

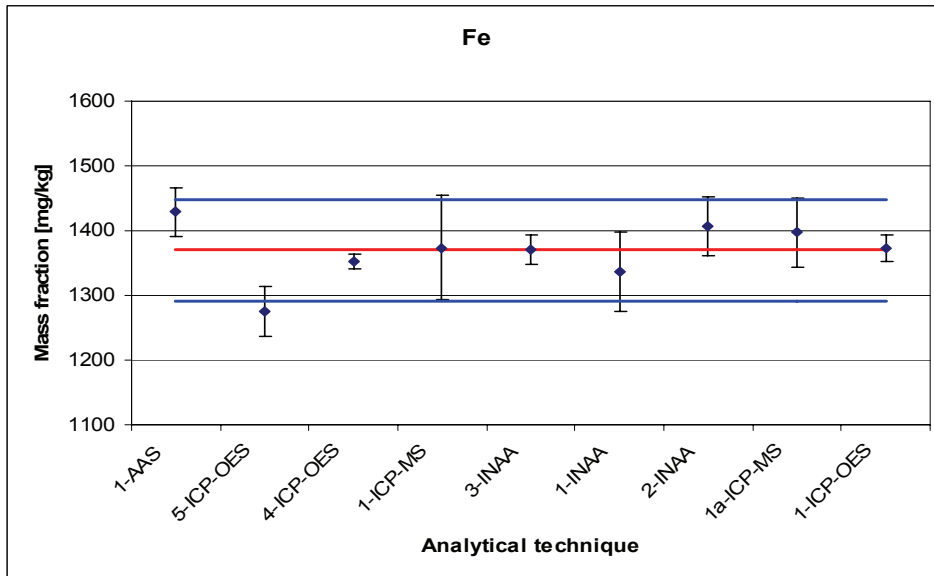


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and information value $\pm U$ for Hg**

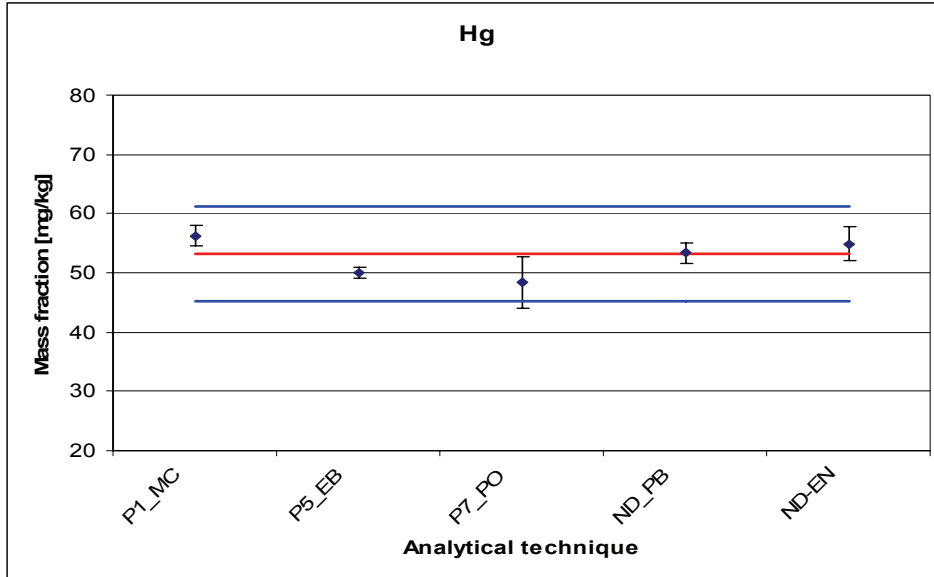


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value $\pm U$ for K**

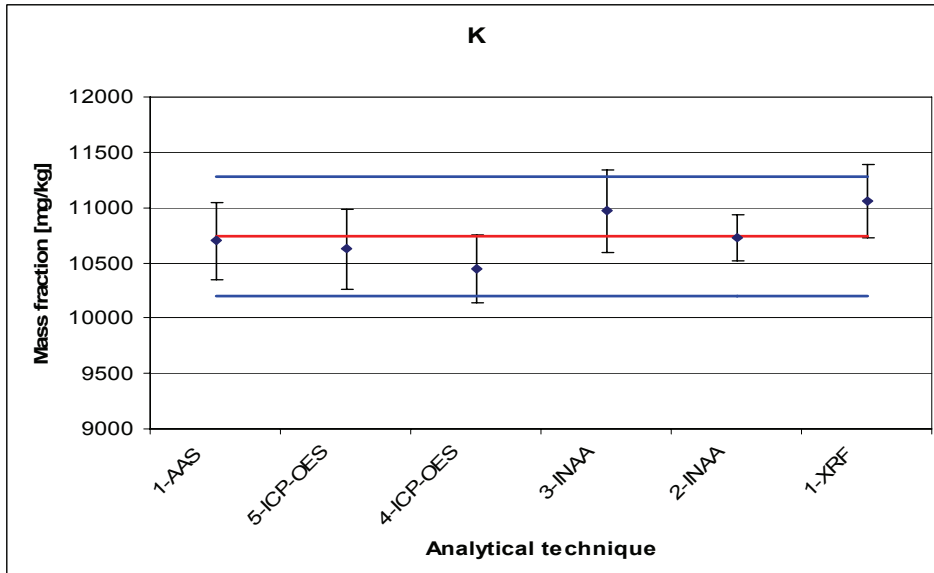


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Mg**

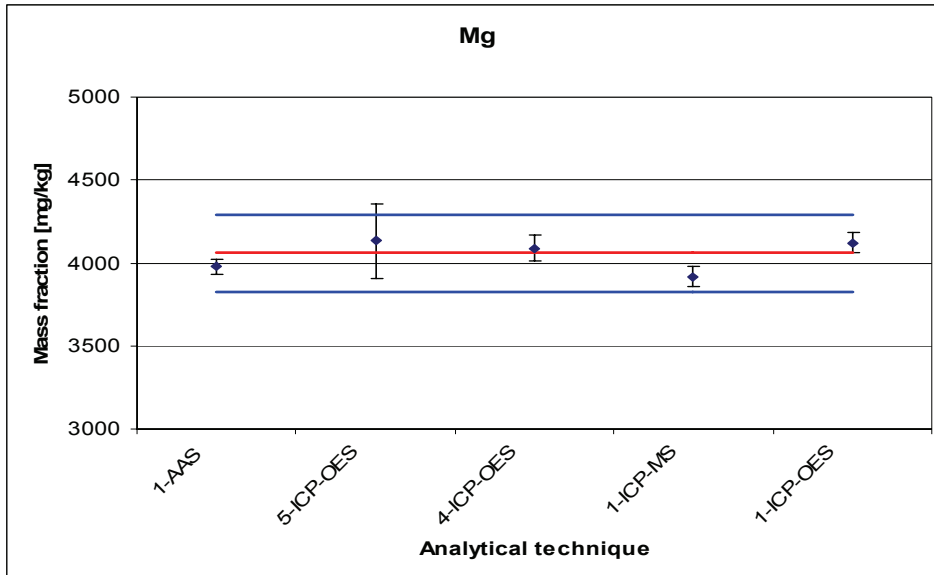


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Mn**

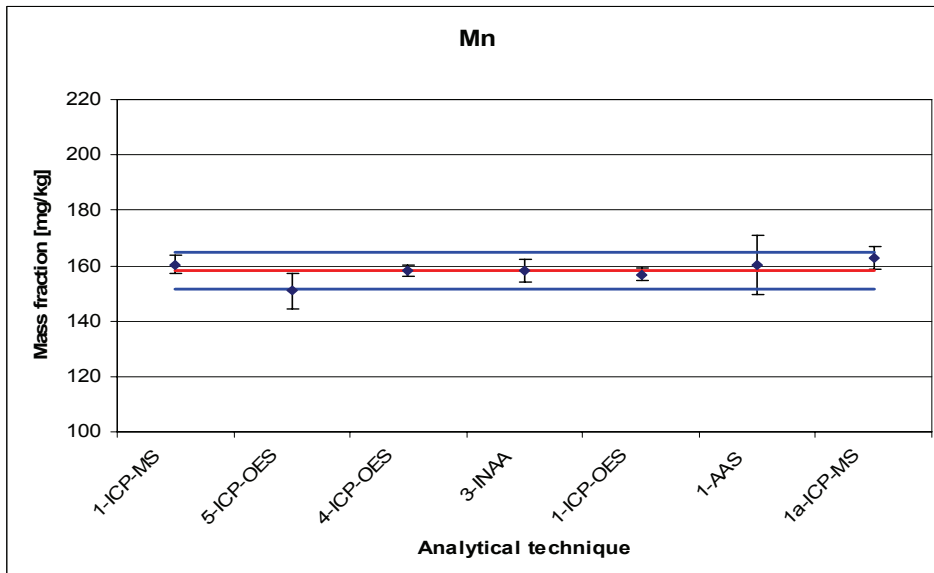


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Na**

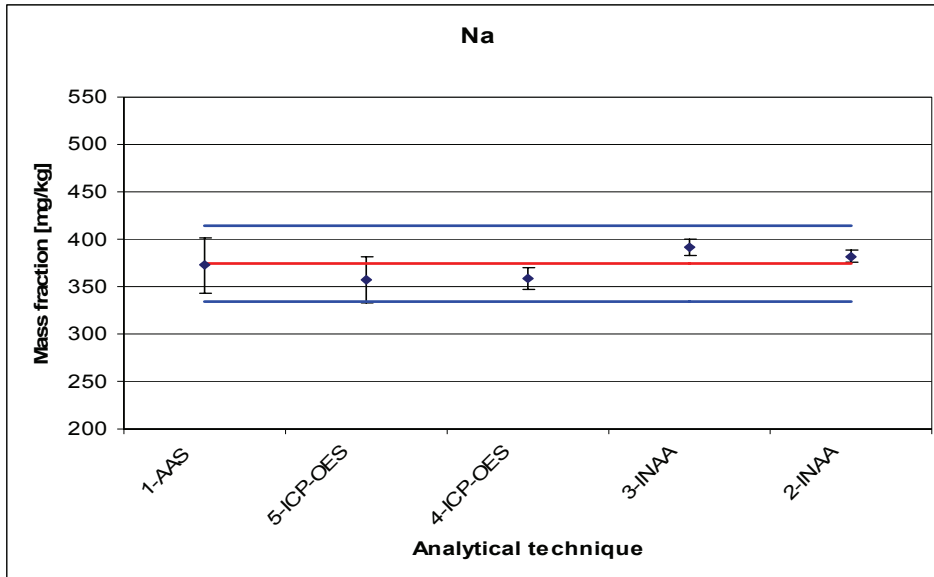


**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Ni**

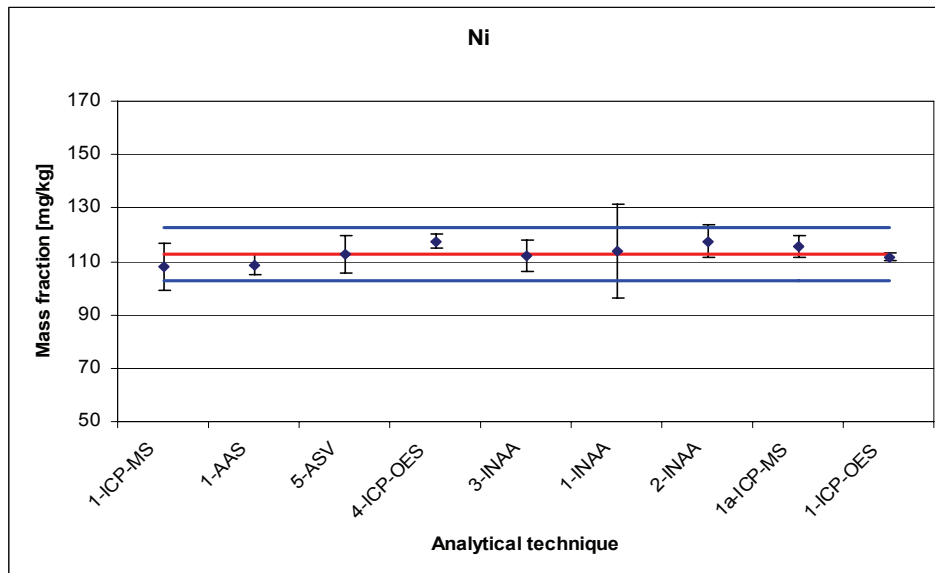

**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Pb**

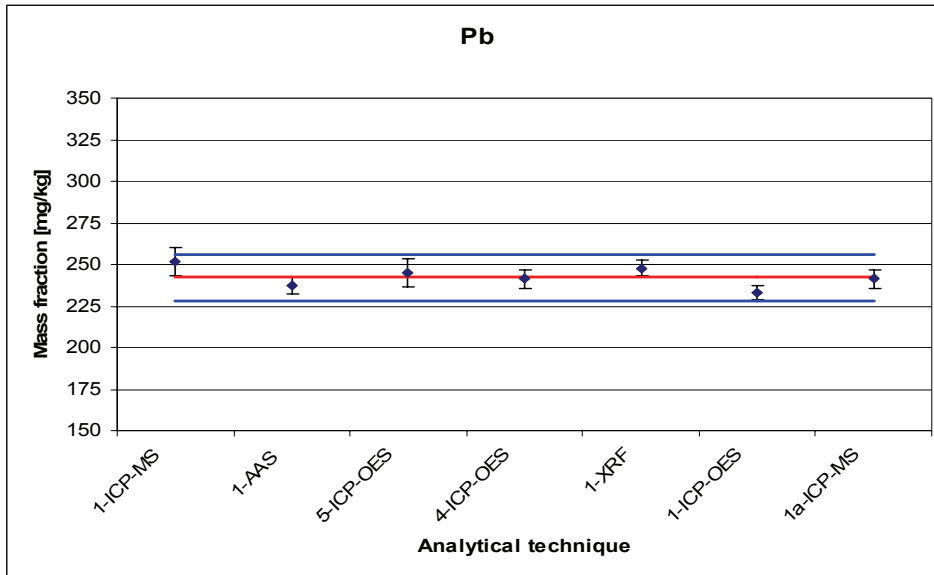

**IAEA 413, Major, minor and trace elements in Algae
Laboratory means and reference value \pm U for Zn**

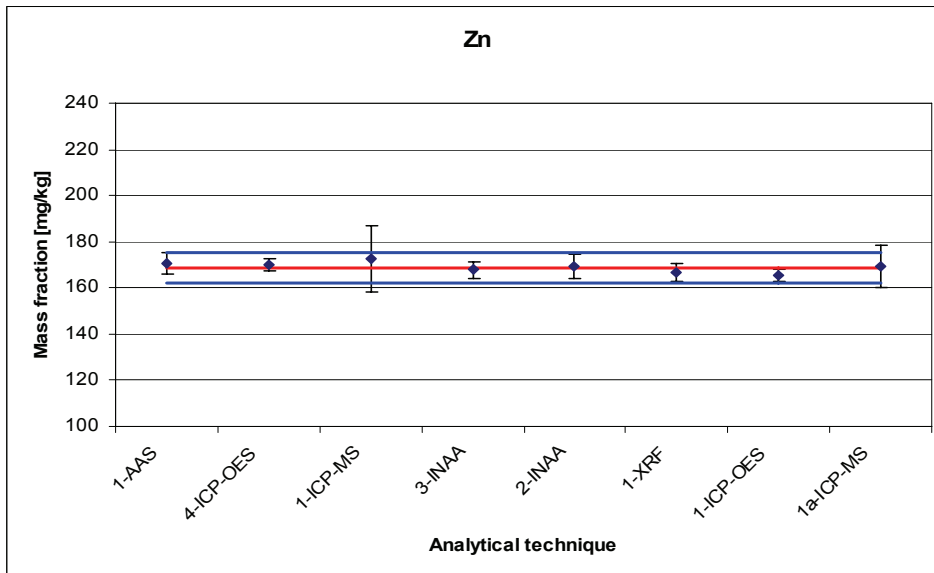

REFERENCES

- [1] INTERNATIONAL STANDARDS ORGANIZATION, ISO Guide 34: General requirements for the competence of reference material producers, Second edition, Geneva, Switzerland, (2000).
- [2] INTERNATIONAL STANDARDS ORGANIZATION, ISO Guide 35: Reference materials - General and statistical principles for certification, Geneva, Switzerland, (2005).
- [3] ZEISLER, R; DEKNER, R, ZEILLER, E, DOUCHA, J, MADER, P, KUČERA, J,: Single cell green algae reference materials with managed levels of heavy metals, *Fresenius J. Anal Chem.* **360** (1996) 429-432.
- [4] DOUCHA, J, LIVANSKY, K, Novel outdoor thin layer high density microalgae culture system: Productivity and operational parameters, *Arch. Hydrobiol.* **106/Algolog. Stud.** **76** (1995) 129-147.
- [5] INTERNATIONAL ATOMIC ENERGY AGENCY, Reference Materials for Microanalytical Nuclear Techniques, IAEA-TECDOC-1295, IAEA, Vienna, Austria, (2002).
- [6] ROSSBACH, M, ZEILLER, E, Assessment of element specific homogeneity in reference materials using microanalytical techniques, *Anal Bioanal Chem.* **377** (2003) 334-339.
- [7] INGAMELLS, C, O, SWITZER P, A proposed sampling constant for use in geochemical analysis, *Talanta* **20** (1973) 547.
- [8] PAUWELS, J, VANDECASTEELE, C, Determination of the minimum sample mass of a solid CRM to be used in chemical analysis, *Fresenius J. Anal. Chem.*, **345** (1993) 121-123.
- [9] STOEPLER, M, KURFUERST, U, GROBECKER, K, H, Der Homogenitaetsfaktor als Kenngrösse fuer pulverisierte Feststoffproben, *Fresenius' J. Anal Chem.*, **322** (1985), 687-691.
- [10] INTERNATIONAL ATOMIC ENERGY AGENCY, Good Analytical Practice in the Analysis of IAEA Materials for Certification, IAEA/AL/096, IAEA, Vienna, Austria, (1996)
- [11] BLAAUW, M, The k_0 -consistent IRI Gamma-Ray catalogue for INAA, Delft, Netherlands, <http://www.tnw.tudelft.nl/> (1996) 141.
- [12] HORWITZ, W, Evaluation of Analytical Methods for Regulation of Foods and Drugs, *Anal. Chem.*, **54** (1982) 67A-76A.
- [13] ZEILLER, E, BENETKA, E, KOLLER, M, SCHORN, R,: Dry mass determination - which role does it play in combined measurement uncertainty? A case study using IAEA 392 and IAEA 413 algae reference materials, *Accred Qual Assur*, **12** (2007) 295-302.
- [14] DUEWER, D, L, A Robust Approach for the Determination of CCQM Key Comparison Values and Uncertainties, Paper presented at the 10th meeting of the CIPM Consultative Committee for Amount of Substance — Metrology in Chemistry, Sevres, France, (2004). www.bipm.info/cc/ccqm/allowed/10/ccqm04-15.pdf.

ABBREVIATIONS

| | |
|---------|--|
| AAS | Atomic absorption spectrometry |
| ASV | Anodic stripping voltammetry |
| CV | Cold vapour |
| ET | Electro thermal |
| ED-XRF | Energy dispersive X-Ray fluorescence |
| HG | Hydride generation |
| HDP | High pressure decomposition |
| IAEA | International Atomic Energy Agency |
| ICP-MS | Inductively coupled plasma mass spectrometry |
| INAA | Instrumental neutron activation analysis |
| ICP-OES | Inductively coupled plasma optical emission spectrometry |
| LDP | Low pressure decomposition |
| ND | Non destructive |
| PMW | Pressurized microwave |
| PTFE | Teflon |
| XRF | X-Ray fluorescence |

SAMPLE PREPARATION CODES

| | |
|------|---|
| ND | Non destructive |
| P1 | PMW, HNO ₃ |
| P1HF | PMW, HNO ₃ , HF |
| P2 | PMW, HNO ₃ , H ₂ O ₂ |
| P3 | PMW, HNO ₃ , HClO ₄ |
| P4 | open, HNO ₃ , HClO ₄ |
| P5 | LPD, HNO ₃ |
| P6 | HPD at 300°C, HNO ₃ |
| P7 | reflux, HNO ₃ |
| P8 | PTFE at 180°C, HNO ₃ |

ABBREVIATIONS IN THE EQUATIONS AND TABLES (number of equation and tables in brackets) if not explained in the text

| | |
|---------------|--|
| b_1 | <i>Slope</i> (Table 3) |
| b_0 | <i>Intercept</i> (Table 3) |
| k | Coverage factor (Table 7) |
| MS_{among} | Mean square (ANOVA) between bottles (2) |
| MS_{within} | Mean square (ANOVA) within bottles (1, 2, 3) |

| | |
|-----------------------|--|
| n | Number of observations (3, Table 3) |
| n_b | Number of bottles (7, Table2) |
| n_i | Number of individual measurement results (7, Table 6) |
| n_L | Number of laboratories participating in the characterization (7, Table 6 and 7) |
| n_0 | (Effective) number of (sub) group members (for complete data sets $n=n_0$) (2, 3) |
| n_R | Number of replicates per bottle (Table2) |
| $RSD\%$ | Relative standard deviation (Table1) |
| s | Standard deviation |
| s_{b1} | Uncertainty of the slope (Table 3) |
| s_{b0} | Uncertainty of the intercept (Table 3) |
| s.e.m | Standard error of the mean (Table 2 and 6) |
| s_i | standard deviation of all individual measurement results (Table 6) |
| s_{bb}^2 | Variance between bottles (2) |
| $S(MM\text{-}median)$ | Standard deviation of Mixture Model Median |
| s_{bL}^2 | Variance between laboratories |
| s_{wb}^2 | Variance within bottles (1) |
| u_x | Combined standard uncertainty of property value (grand mean) (Table 7) |
| u_{wb} | Uncertainty related to within bottle variations (Table 7) |
| u_{bb} | Uncertainty related to between bottle variations (Table 7) |
| u^*_{bb} | Uncertainty related to between bottle heterogeneity (3, Table 2) |
| u_{bL} | Uncertainty related to between laboratory variations (Table 7) |
| U | Expanded uncertainty (coverage factor 2 for 95% probability) (Table 7) |
| $V_{MS_{within}}$ | Degree of freedom of mean square within bottles (3) |
| x_i | Individual result |

PARTICIPATING LABORATORIES

Laboratorio de Radioisotopos, Centro de Energia Nuclear na Agricultura
Universidada de Sao Paulo
Caixa Postal 96
13400-970 Piracicaba
S.P, Brazil
Dr. E. de Nadai Fernandes

KFZ Juelich, Abteilung IPC
Postfach 1913
D-52425 Jülich
Germany
Dr. P. Ostapczuk

Austrian Research Center Seibersdorf
Hauptabteilung Umwelt-Analytik
A-2444 Seibersdorf
Austria
Mag. E. Benetka

Delft University of Technology
Mekelweg 15
NL-2629 JB Delft
The Netherlands
Dr. Peter Bode

International Atomic Energy Agency
Chemistry Unit, Terrestrial Environment Laboratory
IAEA Environment Laboratories
Austria
Mr. M. Campbell, Mr. A. Törvényi, Mr. M. Makarewicz, Ms. E. Zeiller, Mr. S. Bamfort

CONTRIBUTORS TO DRAFTING AND REVIEW

Fajgelj, A. International Atomic Energy Agency

Shakhashiro, A. International Atomic Energy Agency

Törvenyi, A. International Atomic Energy Agency

Zeiller, E. International Atomic Energy Agency



Where to order IAEA publications

In the following countries IAEA publications may be purchased from the sources listed below, or from major local booksellers. Payment may be made in local currency or with UNESCO coupons.

Australia

DA Information Services, 648 Whitehorse Road, Mitcham Victoria 3132
Telephone: +61 3 9210 7777 • Fax: +61 3 9210 7788
Email: service@dadirect.com.au • Web site: <http://www.dadirect.com.au>

Belgium

Jean de Lannoy, avenue du Roi 202, B-1190 Brussels
Telephone: +32 2 538 43 08 • Fax: +32 2 538 08 41
Email: jean.de.lannoy@infoboard.be • Web site: <http://www.jean-de-lannoy.be>

Canada

Bernan Associates, 4611-F Assembly Drive, Lanham, MD 20706-4391, USA
Telephone: 1-800-865-3457 • Fax: 1-800-865-3450
Email: order@bernan.com • Web site: <http://www.bernan.com>

Renouf Publishing Company Ltd., 1-5369 Canotek Rd., Ottawa, Ontario, K1J 9J3
Telephone: +613 745 2665 • Fax: +613 745 7660
Email: order.dept@renoufbooks.com • Web site: <http://www.renoufbooks.com>

China

IAEA Publications in Chinese: China Nuclear Energy Industry Corporation, Translation Section, P.O. Box 2103, Beijing

Czech Republic

Suweco CZ, S.R.O. Klecakova 347, 180 21 Praha 9
Telephone: +420 26603 5364 • Fax: +420 28482 1646
Email: nakup@suweco.cz • Web site: <http://www.suweco.cz>

Finland

Akateeminen Kirjakauppa, PL 128 (Keskuskatu 1), FIN-00101 Helsinki
Telephone: +358 9 121 41 • Fax: +358 9 121 4450
Email: akatilaus@akateeminen.com • Web site: <http://www.akateeminen.com>

France

Form-Edit, 5, rue Janssen, P.O. Box 25, F-75921 Paris Cedex 19
Telephone: +33 1 42 01 49 49 • Fax: +33 1 42 01 90 90 • Email: formedit@formedit.fr

Lavoisier SAS, 14 rue de Provigny, 94236 Cachan Cedex
Telephone: + 33 1 47 40 67 00 • Fax +33 1 47 40 67 02
Email: livres@lavoisier.fr • Web site: <http://www.lavoisier.fr>

Germany

UNO-Verlag, Vertriebs- und Verlags GmbH, August-Bebel-Allee 6, D-53175 Bonn
Telephone: +49 02 28 949 02-0 • Fax: +49 02 28 949 02-22
Email: info@uno-verlag.de • Web site: <http://www.uno-verlag.de>

Hungary

Librotrade Ltd., Book Import, P.O. Box 126, H-1656 Budapest
Telephone: +36 1 257 7777 • Fax: +36 1 257 7472 • Email: books@librotrade.hu

India

Allied Publishers Group, 1st Floor, Dubash House, 15, J. N. Heredia Marg, Ballard Estate, Mumbai 400 001,
Telephone: +91 22 22617926/27 • Fax: +91 22 22617928
Email: alliedpl@vsnl.com • Web site: <http://www.alliedpublishers.com>

Bookwell, 24/4800, Ansari Road, Darya Ganj, New Delhi 110002
Telephone: +91 11 23268786, +91 11 23257264 • Fax: +91 11 23281315
Email: bookwell@vsnl.net • Web site: <http://www.bookwellindia.com>

Italy

Libreria Scientifica Dott. Lucio di Biasio "AEIOU", Via Coronelli 6, I-20146 Milan
Telephone: +39 02 48 95 45 52 or 48 95 45 62 • Fax: +39 02 48 95 45 48

Japan

Maruzen Company, Ltd., 13-6 Nihonbashi, 3 chome, Chuo-ku, Tokyo 103-0027
Telephone: +81 3 3275 8582 • Fax: +81 3 3275 9072
Email: journal@maruzen.co.jp • Web site: <http://www.maruzen.co.jp>

Korea, Republic of

KINS Inc., Information Business Dept. Samho Bldg. 2nd Floor, 275-1 Yang Jae-dong SeoCho-G, Seoul 137-130
Telephone: +02 589 1740 • Fax: +02 589 1746
Email: sj8142@kins.co.kr • Web site: <http://www.kins.co.kr>

Netherlands

Martinus Nijhoff International, Koraalrood 50, P.O. Box 1853, 2700 CZ Zoetermeer
Telephone: +31 793 684 400 • Fax: +31 793 615 698 • Email: info@nijhoff.nl • Web site: <http://www.nijhoff.nl>

Swets and Zeitlinger b.v., P.O. Box 830, 2160 SZ Lisse
Telephone: +31 252 435 111 • Fax: +31 252 415 888 • Email: info@swets.nl • Web site: <http://www.swets.nl>

New Zealand

DA Information Services, 648 Whitehorse Road, MITCHAM 3132, Australia
Telephone: +61 3 9210 7777 • Fax: +61 3 9210 7788
Email: service@dadirect.com.au • Web site: <http://www.dadirect.com.au>

Slovenia

Cankarjeva Založba d.d., Kopitarjeva 2, SI-1512 Ljubljana
Telephone: +386 1 432 31 44 • Fax: +386 1 230 14 35
Email: import.books@cankarjeva-z.si • Web site: <http://www.cankarjeva-z.si/uvoz>

Spain

Díaz de Santos, S.A., c/ Juan Bravo, 3A, E-28006 Madrid
Telephone: +34 91 781 94 80 • Fax: +34 91 575 55 63 • Email: compras@diazdesantos.es
carmela@diazdesantos.es • barcelona@diazdesantos.es • julio@diazdesantos.es
Web site: <http://www.diazdesantos.es>

United Kingdom

The Stationery Office Ltd, International Sales Agency, PO Box 29, Norwich, NR3 1 GN
Telephone (orders): +44 870 600 5552 • (enquiries): +44 207 873 8372 • Fax: +44 207 873 8203
Email (orders): book.orders@tso.co.uk • (enquiries): book.enquiries@tso.co.uk • Web site: <http://www.tso.co.uk>

On-line orders:

DELTA Int. Book Wholesalers Ltd., 39 Alexandra Road, Addlestone, Surrey, KT15 2PQ
Email: info@profbooks.com • Web site: <http://www.profbooks.com>

Books on the Environment:

Earthprint Ltd., P.O. Box 119, Stevenage SG1 4TP
Telephone: +44 1438748111 • Fax: +44 1438748844
Email: orders@earthprint.com • Web site: <http://www.earthprint.com>

United Nations (UN)

Dept. 1004, Room DC2-0853, First Avenue at 46th Street, New York, N.Y. 10017, USA
Telephone: +800 253-9646 or +212 963-8302 • Fax: +212 963-3489
Email: publications@un.org • Web site: <http://www.un.org>

United States of America

Bernan Associates, 4611-F Assembly Drive, Lanham, MD 20706-4391
Telephone: 1-800-865-3457 • Fax: 1-800-865-3450
Email: order@bernan.com • Web site: <http://www.bernan.com>

Renouf Publishing Company Ltd., 812 Proctor Ave., Ogdensburg, NY, 13669
Telephone: +888 551 7470 (toll-free) • Fax: +888 568 8546 (toll-free)
Email: order.dept@renoufbooks.com • Web site: <http://www.renoufbooks.com>

Orders and requests for information may also be addressed directly to:

Sales and Promotion Unit, International Atomic Energy Agency

Vienna International Centre, PO Box 100, 1400 Vienna, Austria
Telephone: +43 1 2600 22529 (or 22530) • Fax: +43 1 2600 29302
Email: sales.publications@iaea.org • Web site: <http://www.iaea.org/books>

INTERNATIONAL ATOMIC ENERGY AGENCY
VIENNA
ISSN 2074-7659